

Title:

Magnetically controlled capsule endoscopy (MCCE) as an alternative for gastric evaluation

Authors:

Gerardo Blanco-Velasco, Enrique Pérez-Cuadrado Martínez, Pilar Esteban-Delgado, Asunción María Gonzálvez-Gasch

DOI: 10.17235/reed.2025.11597/2025 Link: <u>PubMed (Epub ahead of print)</u>

Please cite this article as:

Blanco-Velasco Gerardo, Pérez-Cuadrado Martínez Enrique, Esteban-Delgado Pilar, Gonzálvez-Gasch Asunción María. Magnetically controlled capsule endoscopy (MCCE) as an alternative for gastric evaluation. Rev Esp Enferm Dig 2025. doi: 10.17235/reed.2025.11597/2025.

This is a PDF file of an unedited manuscript that has been accepted for publication. As a service to our customers we are providing this early version of the manuscript. The manuscript will undergo copyediting, typesetting, and review of the resulting proof before it is published in its final form. Please note that during the production process errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.



Magnetically controlled capsule endoscopy (MCCE) as an alternative for gastric evaluation

Gerardo Blanco-Velasco¹, Enrique Francisco Pérez-Cuadrado-Martínez², Pilar Esteban-Delgado², Asunción María Gonzálvez-Gasch³.

¹Digestive Diseases Service, Hospital Quironsalud Torrevieja, Torrevieja, Alicante

²Digestive Diseases Service, Hospital General Universitario Morales Meseguer, Murcia, Murcia

³Internal Medicine Service, Hospital Quironsalud Torrevieja, Torrevieja, Alicante

Corresponding Author:

Gerardo Blanco-Velasco

E-mail: gerardoblancov@hotmail.com

Address: Calle Almirante Malaspina 3, 2C, Murcia, Murcia, España, CP 30007.



Conventional capsule endoscopy (CE) is widely used for examining the gastrointestinal tract, especially the small bowel (SB) and colon.[1] Magnetically Controlled Capsule Endoscopy (MCCE) is an innovative approach for gastric evaluation, suitable for selected cases. The system utilizes a magnetically controlled capsule operated by an external magnetic robotic C-arm, with two joysticks providing precise control orientation. For optimal stomach visualization, the patient drinks at least 1 liter of water to float the capsule and distend the stomach, replacing air insufflation. The operator can then guide the MCCE for a comprehensive 360-degree assessment, aided by adjusting the patient's position, or use automatic mode, where the patient simply lies down after swallowing. In both modes, artificial intelligence detects key gastric landmarks (cardia, fundus, incisura, antrum, and pylorus).[2,3] MCCE offers detailed examination of both the esophagus and SB. The external magnetic control can stop the capsule in the esophagus after swallowing, facilitating the inspection of the esophageal mucosa. After visualizing the esophagus and stomach, the capsule naturally moves into the SB, functioning similarly to standard CE examining the entire SB.

We included 30 patients requiring SB evaluation (bleeding, Crohn's disease, anemia, diarrhea) in the study, excluding those with pacemakers or metal prostheses. All procedures were performed at Quironsalud Torrevieja Hospital using the NaviCam MCCE® (Ankon Technologies, Shanghai) from January to August 2025. One patient, who was unable to swallow the capsule, had it placed endoscopically. For the other 29 patients, the NaviCam Xpress robotic arm enabled comprehensive visualization of the esophagus and stomach. The median patient age was 63 years (IQR: 48-70). Median transit times were 83 seconds (IQR 20-155) for the esophagus, 50 minutes (IQR 34-75) for the stomach, and 305 minutes (IQR 233-433) for the SB. Complete, manually guided gastroscopy with MCCE took a median of 19 minutes (IQR 18-22). In eight cases (27.5%), the MCCE was advanced into the duodenum using external magnetic control. Gastric findings are detailed in figure 1, with additional diagnoses of esophagitis in two patients and a hiatal hernia in another. (Fig. 2) MCCE revealed previously unknown findings in twelve patients, changing the management in two (gastric ulcers and esophagitis).



Although gastroscopy remains the gold standard, MCCE's could be considered as an alternative in certain indications, as contraindication or refusal to standard gastroscopy, screening (cancer, varices, etc.), gastroscopies without biopsy, or pandemic situations must be evaluated with further investigations. To our knowledge, this is the first study using MCCE in Spain.

References

- 1.- Carretero C, Prieto de Frías C, Angós R, Betés M, Herráiz M, et al. Pan-enteric capsule for bleeding high-risk patients. Can we limit endoscopies? Rev Esp Enferm Dig. 2021;113:580-584.
- 2.- Zhang Y, Zhang Y, Huang X. Development and Application of Magnetically Controlled Capsule Endoscopy in Detecting Gastric Lesions. Gastroenterol Res Pract 2021;2021:2716559.
- 3.- Szalai M, Helle K, Lovász BD, Finta Á, Rosztóczy A, et al. First prospective European study for the feasibility and safety of magnetically controlled capsule endoscopy in gastric mucosal abnormalities. World J Gastroenterol 2022;28:2227-2242.

Figure 1. Outcomes of magnetically controlled capsule endoscopy in the stomach

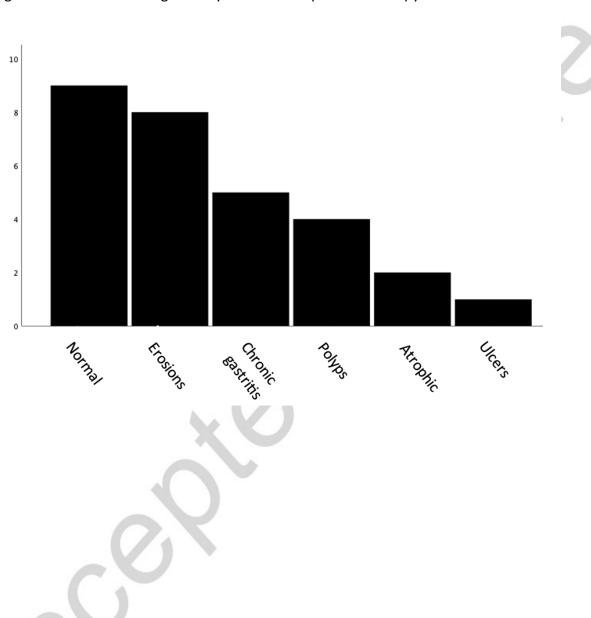




Figure 2. Collage of representative images obtained with Magnetically Controlled Capsule Endoscopy. Normal esophageal landmarks: a) upper part, b) lower part, c) Z-line; and pathological esophageal finding: d) esophagitis). Normal gastric landmarks: e) fundus, f-g) body, h) incisura, i) antrum, j) pylorus); and gastric pathology: k) atrophic gastritis, l) polyp. Small bowel normal images: m) ampulla of Vater (blue arrow), n) jejunum; and small bowel pathology: o) angioectasia, p) ulcer.

